



AG-375M

Heat & Massage



Whiskey




Java



Otter




Rose

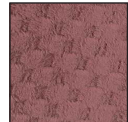

Champagne



Ebony Vinyl

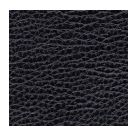

shown in Whiskey

This lift and recline chair comes standard with an overstuffed seam back and our heater/vibrating massage unit. Featuring a full chaise pad for comfort and leg support, the 375M-HM has a 375 lb. weight capacity and is ideal for someone 5'4" to 5'10".



AG-375M-HM

Model Number	375M-HM
Size	Medium
Weight Capacity	375 LB
Recline Positions	3
Standard Fabrics	5
Chaise Seat	Yes
Back Type	Seam
Battery Back Up	Yes
Overall Width	33"
Overall Height	41"
Width Between Arms	21"
Floor to Top of Seat	19"
Seat Depth	21"
Seat to Top of Back	28"
Seat to Top of Arm	7.5"
Floor to Back Edge of Seat Raised	27"
Floor to Front Edge of Seat Raised	24"
Seat Angle Raised	15°
Floor to Top of Foot Rest Reclined	22"
Floor to Top of Back Reclined	31"
Back Angle Closed	120°
Back Angle Reclined	155°
Foot Rest Angle Closed	90°
Foot Rest Angle Opened	170°
Distance Required from Wall Reclined	17"
Extension Reclined from Back of Seat to End of Foot Rest	42"
Lifetime Warranty, Chair Frame, Lift Frame, Recline Mechanism	Yes
Warranty on Electrical parts	1 YR
Warranty on Mechanical Labor	No
Warranty on Service	No

* See Customize Chart for dimension limits
 Sizes may vary depending on fabric, filling material, upholstery and carpet thickness.
 All measurements taken on concrete floor using levels and metal rulers.

User-friendly feature of the AG-375M-HM lift chair:

The 375M-HM recliner has a battery backup system in the electrical transformer for two 9 volt batteries. In the event of a power failure, the battery system will raise the chair from the reclined to the seated position. The batteries should be replaced after one use.



The AmeriGlide 375M-HM is constructed of the highest quality materials. The frame is composed of northern maple hardwood, poplar and engineered furniture laminate products. The chair's frame joints are joined using staples, glue and mortise and tenon joinery for long lasting durability.

The sinuous spring system on the AmeriGlide 375M-HM follows the manufacturing practices found in the nation's finest shops since its introduction in 1934. Traditional furniture building techniques are the foundation of the AmeriGlide 375M-HM. The springs are secured in the front and back with insulated clips to prevent squeaking and sagging. Our seats have four seat springs. To further enhance the strength and durability of our springs, we secure them with tie wire and material to the seat side rail. We follow the same process with the springs found in the backrest, further reducing noise and sagging for a longer lifetime of use. This forms the foundation for attaching 1.8 density seat foam with California 117 firerating.

AmeriGlide is proud to be a pioneer in developing a chaise pad that is superior in the marketplace. The pad is sewn to a seat covering and then secured to the ottoman, preventing movement of the chaise pad while the chair is in use. The chaise pad is further secured from the front to the back and then side to side. The chaise pad is then apex sewed, which further secures the fabric and foam to the seat for years of trouble-free use.



AmeriGlide
 Accessibility
 Solutions



Literature is current at time of printing. AmeriGlide reserves the right to make changes to the product or literature at any time.



Complies with UL 73 and CSA C22.2 No. 68

Listing # E112417

